GE Healthcare

# Power Upgrade Option Installation Manual

OPERATING DOCUMENTATION

5428387-1EN

Rev 2

© 2014 General Electric Company.

All rights reserved.

---

# **LEGAL NOTES**

---

## **TRADEMARKS**

Adobe, the Adobe logo, Acrobat, the Acrobat logo, Exchange, and PostScript are trademarks of Adobe Systems Incorporated or its subsidiaries and may be registered in certain jurisdictions.

Microsoft is a registered trademark and Windows is a trademark of Microsoft Corporation.

All other products and their name brands are trademarks of their respective holders.

## **COPYRIGHTS**

All Material Copyrighted (c) 2014 by the GE Company, All rights reserved.

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## **DAMAGE IN TRANSPORTATION**

All packages should be closely examined at time of delivery. If damage is apparent, have notation "damage in shipment" written on all copies of the freight or express bill before delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 8.
- Fill out a report on <http://us44hdd21/sctq/InstallFulfill/InstalFulfillment.htm>
- Contact your local service coordinator for more information on this process.

Rev. Jan. 5, 2005

## **CERTIFIED ELECTRICAL CONTRACTOR STATEMENT**

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

## **IMPORTANT...X-RAY PROTECTION**

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Medical Systems Group, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protection, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Medical Systems Group, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

Various protective materials and devices are available. It is urged that such materials or devices be used.

## LITHIUM BATTERY CAUTIONARY STATEMENTS

### CAUTION

**Risk of Explosion.**

**Danger of explosion if battery is incorrectly replaced.**

**Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer's instructions.**

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## OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the Healthcare PQR Process to report all omissions, errors, and defects in this publication.

---

# Revision History

---

Revision	Date	Reason for change
2	May 22, 2014	Added Information of Revolution EVO
1	October 26, 2011	Initial Release

---

# List of Effected Pages

---

PAGES	REVISION	PAGES	REVISION	PAGES	REVISION
1 through 36	2				

# Table of Contents

## Preface

**Publication Conventions ..... 19**

**Section 1.0**

**Safety & Hazard Information ..... 19**

1.1 Text and Character Representation..... 19  
1.2 Graphical Representation ..... 20

**Section 2.0**

**Publication Conventions ..... 21**

2.1 General Paragraph and Character Styles..... 21  
2.2 Page Layout..... 21  
2.3 Computer Screen Output/Input Text Character Styles ..... 22  
2.4 Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)..... 22

## Chapter 1

**Introduction..... 23**

## Chapter 2

**Power upgrade from 33KW to 48KW (50KVA > 75kVA) ..... 25**

**Section 1.0**

**Personnel Requirements ..... 25**

**Section 2.0**

**Overview ..... 25**

2.1 Flowchart ..... 25

**Section 3.0**

**Preliminary Requirements..... 26**

3.1 Tools and Test Equipment..... 26  
3.2 Consumables..... 26  
3.3 Material List ..... 26  
3.4 Safety..... 26  
3.5 Required Conditions ..... 26

**Section 4.0**

**Installation Procedure..... 27**

4.1 Remove 50KVA Power Option Name Plates and MHLW option plate..... 27  
4.2 Attach 75KVA Power Option Name Plates ..... 28  
4.3 Attach MHLW Option Plate..... 28  
4.4 Install 75kVA Option SW Key ..... 29

## Chapter 3

<b>Power upgrade from 48KW to 72KW (75KVA &gt; 100kVA)</b> .....	<b>31</b>
<b>Section 1.0</b>	
<b>Personnel Requirements</b> .....	<b>31</b>
<b>Section 2.0</b>	
<b>Overview</b> .....	<b>31</b>
2.1 Flowchart.....	31
<b>Section 3.0</b>	
<b>Preliminary Requirements</b> .....	<b>32</b>
3.1 Tools and Test Equipment .....	32
3.2 Consumables .....	32
3.3 Material List.....	32
3.4 Safety .....	32
3.5 Required Conditions.....	32
<b>Section 4.0</b>	
<b>Installation Procedure</b> .....	<b>33</b>
4.1 Remove 75kVA Power Option Name Plates .....	33
4.2 Attach 100kVA MHLW Option Plate.....	34
4.3 Install 100kVA Option SW Key.....	34

---

# Preface

## Publication Conventions

---

Purpose: Standardize conventions for representing information in a uniform way of communicating information to a reader in a consistent manner. Conventions are used so that the reader can easily recognize the actions or decisions that must be made. There are a number of character and paragraph styles used in this publication to accomplish this task. Please become familiar with them before proceeding forward.

## Section 1.0

### Safety & Hazard Information

#### 1.1 Text and Character Representation

Within this publication, different paragraph and character styles have been used to indicate potential hazards. Paragraph prefixes, such as hazard, caution, danger and warning, are used to identify important safety information. Text (Hazard) styles are applied to the paragraph contents that is applicable to each specific safety statement. Words describe the type of potential hazard that may be encountered and are placed immediately before the paragraph it modifies. Safety information will normally include:

- Type of potential Hazard
- Nature of potential injury
- Causative condition
- How to avoid or correct the causative condition

#### EXAMPLES OF HAZARD STATEMENTS USED

It's important that you read and understand hazard statements, and not just ignore them.

**DANGER  
SEVERE  
PERSONAL  
INJURY**

**DANGER IS USED WHEN A HAZARD EXISTS WHICH WILL CAUSE SEVERE PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:**

- **ELECTROCUTION**
- **CRUSHING**
- **RADIATION**

**WARNING  
SERIOUS  
PERSONAL  
INJURY**

**WARNING IS USED WHEN A HAZARD EXISTS WHICH COULD OR CAN CAUSE SERIOUS PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:**

- **POTENTIAL FOR SHOCK**
- **EXPOSED WIRES**
- **FAILURE TO TAG AND LOCKOUT SYSTEM POWER COULD ALLOW FOR UN-COMMANDED MOTION.**

 **CAUTION**  
**Minor Personal Injury**  
Caution is used when a hazard exists which can or could cause minor injury to self or others if instructions are ignored. They include for example:

- Loss of critical patient data
- Crush or pinch points
- Sharp objects

 **NOTICE**  
**Property Damage Only**  
Notice is used when a hazard is present that can cause property damage but has absolutely no personal injury risk. They can include:

- Disk drive will be destroyed
- Internal mechanical damage, such as to the x-ray tube
- Coasting the rotor through resonance.

## 1.2 Graphical Representation

Important information will always be preceded by the exclamation point  contained within a triangle, as seen throughout the Preface chapter. In addition to text, several different graphical icons (symbols) may be used to make you aware of specific types of hazards that could possibly cause harm.

<p><b>ELECTRICAL</b></p>  <p><b>LASER</b></p>  <p><b>LASER LIGHT</b></p>	<p><b>HAZARDS MECHANICAL</b></p>  <p><b>HEAT</b></p> 	<p><b>RADIATION</b></p>  <p><b>PINCH</b></p> 
<p><b>AVOID STATIC ELECTRICITY</b></p> 	<p><b>PROCEDURAL TAG AND LOCK OUT</b></p> 	<p><b>WEAR EYE PROTECTION</b></p>  <p><b>EYE PROTECTION</b></p>

# Section 2.0 Publication Conventions

## 2.1 General Paragraph and Character Styles

Prefixes are used to highlight important non-safety related information. Paragraph prefixes (such as purpose, example, comment or note) are used to identify important but non-safety related information. Text styles are also applied to text within each paragraph modified by the specific prefix.

### EXAMPLES OF PREFIXES USED FOR GENERAL INFORMATION


**Purpose:** Introduces and provides meaning as to the information contained within the chapter, section or subsection, Such as used at the beginning this chapter for example.

**Note:** Conveys information that should be considered important to the reader.

**Example:** Used to make the reader aware that the paragraph(s) that follow are examples of information possibly stated previously.

**Comment:** Represents "additional" information that may or may not be relevant to your situation.

## 2.2 Page Layout

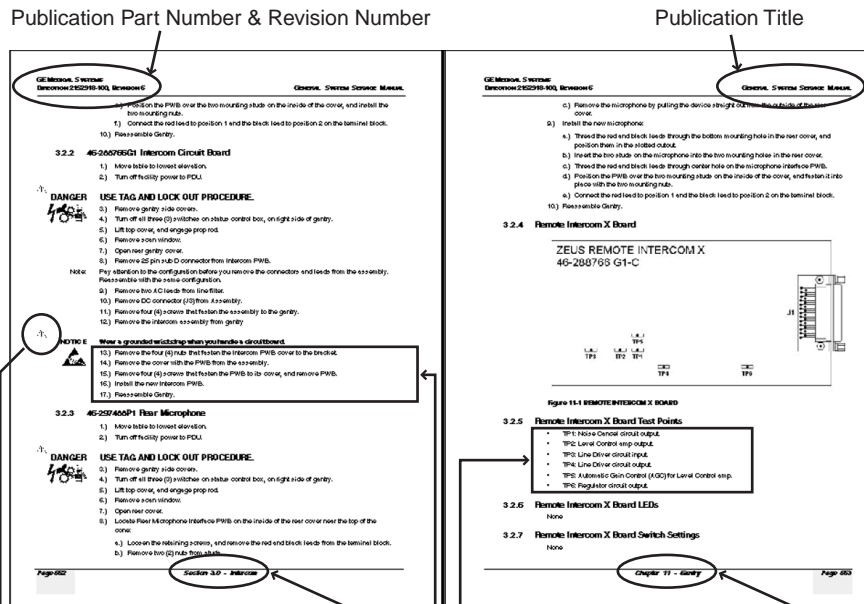The current section and its title are always shown in the footer of the left (even) page.

An exclamation point in a triangle is used to indicate important information to the user.

Paragraphs preceded by **Alphanumeric** characters (e.g. numbers) contain information that must be followed in a **specific order**.

The current chapter and its title are always shown in the footer of the right (odd) page.

Paragraphs preceded by a **symbol** (e.g. bullets) contain information that has **no specific order**.

Headers and footers in this publication are designed to allow you to quickly identify your location. The document's part number and revision number appears in every header on every page. Odd

numbered page footers indicate the current chapter, its title and current page number. Even page footers show the current section and its title, as well current page number.

## 2.3 Computer Screen Output/Input Text Character Styles

Within this publication, mono-spaced character styles (fonts) are used to indicate computer text that's either screen input and output. Mono spaced fonts, such as courier, are used to indicate text direction. When you type at your keyboard, you are generating computer input. Occasionally you will see the math operators "greater than" and "less than" symbols used to indicate the start and finish of variable output. When reading text generated by the computer, you are reading it as computer generated output. In addition to direction, characters are italicized (e.g. *italics*) to indicate information specific to your system or site.

Example:  
Fixed Output

This paragraph's fonts represents computer generated screen "fixed" output. It's output is fixed from the sense that it does not vary from application to application. It's the most commonly used style used to indicate filenames, paths and text that does not change from system to system. The character style used is a fixed width such as courier.

Example:  
Variable Output

*This paragraph's font represents computer screen output that is "variable". Its used to represent output that varies from application to application or system to system. Variable output is sometimes found placed between greater than and lesser than operators for clarification. For example: <variable\_ouput> or <3.45.120.3>. In both cases, the < and > operators are not part of the actual input.*

Example:  
Fixed Input

**This paragraph's font represents fixed input. It's computer input that is typed-in via the keyboard. Typed input that does not vary from application to application or system to system. Fixed text the user is required to supply as input. For example: cd /usr/3p**

Example:  
Variable Input

*This paragraph's font represents computer input that can vary from application to application or system to system. With variable text, the user is required to supply system dependent input or information. Variable input sometimes is placed between greater-than and less-than operators. For example: <variable\_input>. In these cases, the (<>) operators would be dropped prior to input. For example: ypcat hosts | grep <3.45.120.3> would be typed into the computer as ypcat hosts | grep 3.45.120.3 without the greater-than and less-than operators.*

## 2.4 Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)

Different character styles are used to indicate actions requiring the reader to press either a hard or soft button, switch or key. Physical hardware, such as buttons and switches, are called hard keys because they are hard wired or mechanical in nature. A keyboard or on/off switch would be a hard key. Software or computer generated buttons are called soft keys because they are software generated. Software driven menu buttons are an example of such keys. Soft and hard keys are represented differently in this publication.

Example:  
Hard Keys

A power switch **ON/OFF** or a keyboard key like **ENTER** is indicated by applying a character style that uses both over and under-lined bold text that is bold. This is a hard key.

Example:  
Soft Keys

Whereas the computer MENU button that you would click with your mouse or touch with your hand uses over and under-lined regular text. This is a soft key.

**END OF SECTION**

---

# Chapter 1

## Introduction

---

### APPLICABLE PRODUCT

This manual applies to the following products.

- Optima CT660
- Revolution EVO

This manual contains the following upgrade procedures.

- Power Upgrade from 33kW to 48kW (50kVA -> 75kVA) ... Chapter 2
- Power Upgrade from 48kW to 72kW (75kVA -> 100kVA) ... Chapter 3

*This Page Intentionally Blank.*

# Chapter 2

## Power upgrade from 33KW to 48KW (50KVA > 75kVA)


### Section 1.0 Personnel Requirements

Required Persons	Preliminary Reqs	Procedure	Finalization
1	Timing Not Applicable	1 hours	0 mins

### Section 2.0 Overview

This document covers the 48KW (75kVA) power upgrade Installation procedure for 33KW (50kVA) configuration system.

#### 2.1 Flowchart
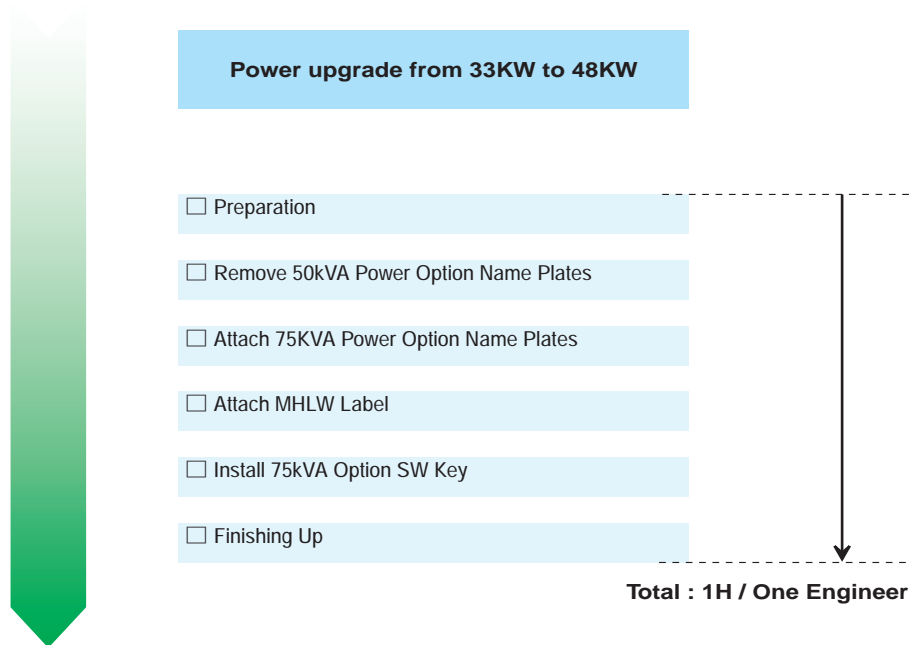
## Section 3.0 Preliminary Requirements

### 3.1 Tools and Test Equipment

None

### 3.2 Consumables

None

### 3.3 Material List

Mark	Number	Qty.	Description
1	5418764	1.0	48kW_75kVA Option SW Key with DVD media
2	5428443	1.0	MHLW Option Plate for Power UPG from 33KW to 48KW
3	5370997	1.0	75kVA Power Option Name Plate on Gantry
4	5371434	1.0	75kVA Power Option Name Plate on PDU
5	5428387-1EN	1.0	Power upgrade Installation Manual (from 33KW to 48KW)

Table 2-1 Power upgrade to 48KW (5418752)

### 3.4 Safety

No general safety information.

### 3.5 Required Conditions

This kit must be installed by GE filled engineer.

## Section 4.0 Installation Procedure

### 4.1 Remove 50kVA Power Option Name Plates and MHLW option plate

- 1.) Remove the 50kVA Power Option Name Plate near the rating plate on the left side of Gantry. Refer to [Figure 2-1](#) ①.
- 2.) Remove the MHLW option plate. Refer to [Figure 2-1](#) ②.

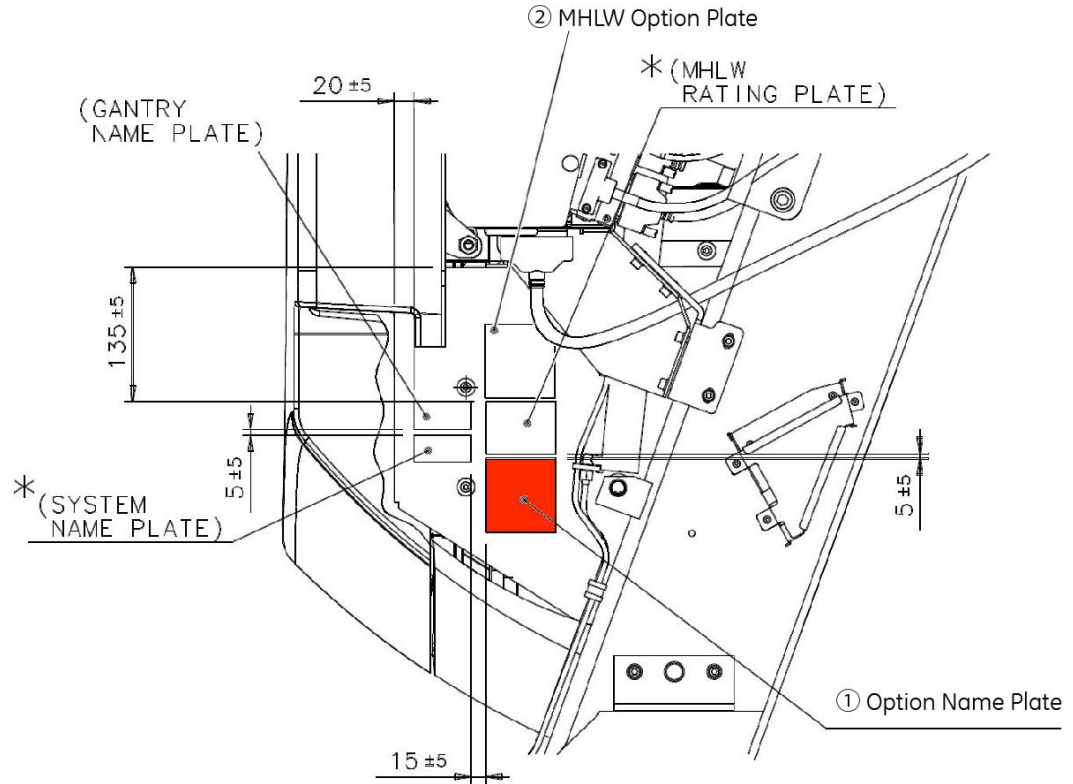
Figure 2-1 Location of Option Name Plate

- 3.) Remove the 50kVA Power Option Name Plate on the PDU rear cover. Refer to [Figure 2-2](#).

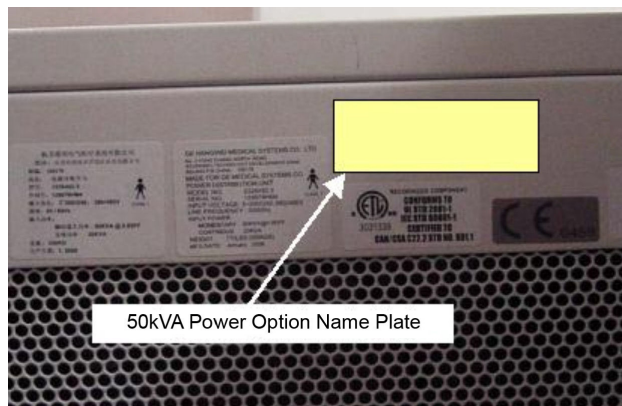
Figure 2-2 Location of Option Name Plate

## 4.2 Attach 75KVA Power Option Name Plates

- 1.) Attach the Gantry 75kVA Power Option Name Plate to the same location where the 50kVA Power Option Name Plate was attached. Refer to [Figure 2-1](#) ①.

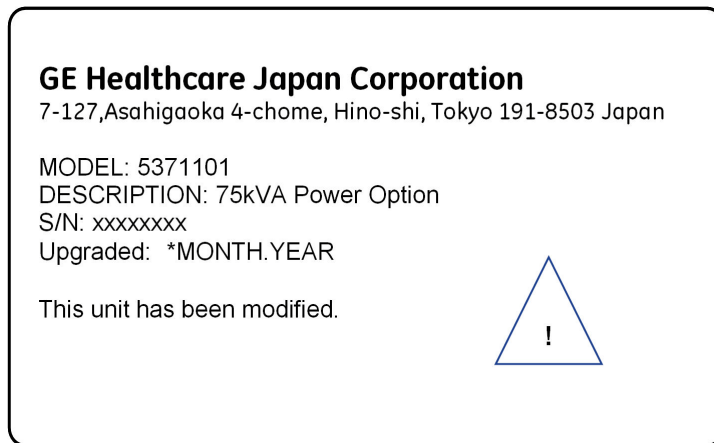
Figure 2-3 75kVA Option Name Plate - Gantry

- 2.) Attach the PDU 75kVA Power Option Name Plate to the same location where the 50kVA Power Option Name Plate was attached. Refer to [Figure 2-2](#).

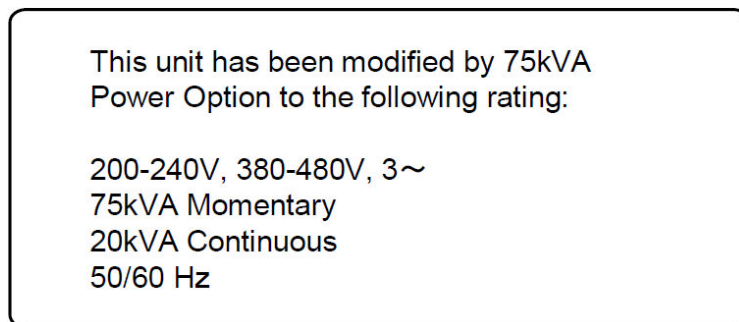
Figure 2-4 75kVA Option Name Plate - PDU

## 4.3 Attach MHLW Option Plate

- 1.) Attach the MHLW Option Plate near the rating plate on the right side of Gantry. Refer to [Figure 2-1](#) ②.

Figure 2-5 MHLW Option Plate

## 4.4 Install 75kVA Option SW Key

- 1.) Install the provided option software key.  
Go to CSD desktop -> select Configuration -> Install Options
- 2.) Reboot the system.

*This Page Intentionally Blank.*

# Chapter 3

## Power upgrade from 48KW to 72KW (75KVA > 100kVA)


### Section 1.0 Personnel Requirements

Required Persons	Preliminary Reqs	Procedure	Finalization
1	Timing Not Applicable	1 hours	0 mins

### Section 2.0 Overview

This document covers the 72KW (100kVA) power upgrade Installation procedure for 48KW (75kVA) configuration system.

#### 2.1 Flowchart
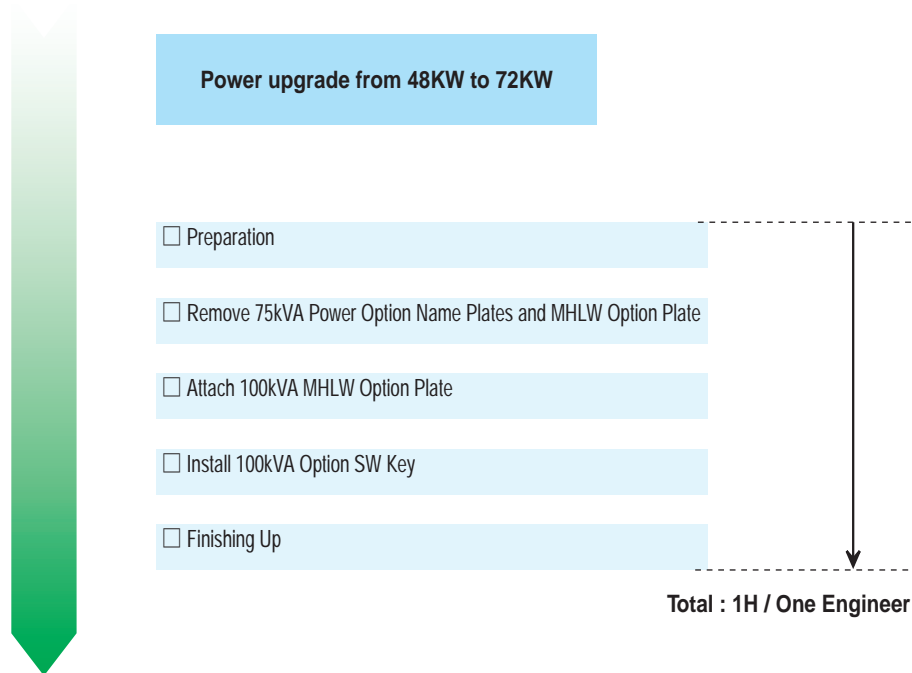
## Section 3.0 Preliminary Requirements

### 3.1 Tools and Test Equipment

None

### 3.2 Consumables

None

### 3.3 Material List

Mark	Number	Qty.	Description
1	5371520	1.0	72kW_100kVA Option SW Key - with DVD
2	5428444	1.0	MHLW Option Plate for Power UPG from 48KW to 72KW
	5428387-1EN	1.0	Power upgrade Installation Manual

**Table 3-1 Power upgrade from 48KW to 72KW (5400970-2)**

### 3.4 Safety

No general safety information.

### 3.5 Required Conditions

This kit must be installed by GE filled engineer.

# Section 4.0 Installation Procedure

## 4.1 Remove 75kVA Power Option Name Plates

- 1.) Remove the 75kVA Power Option Name Plate ① and MHLW Option Plate ② on the left side of Gantry. Refer to [Figure 3-1](#).

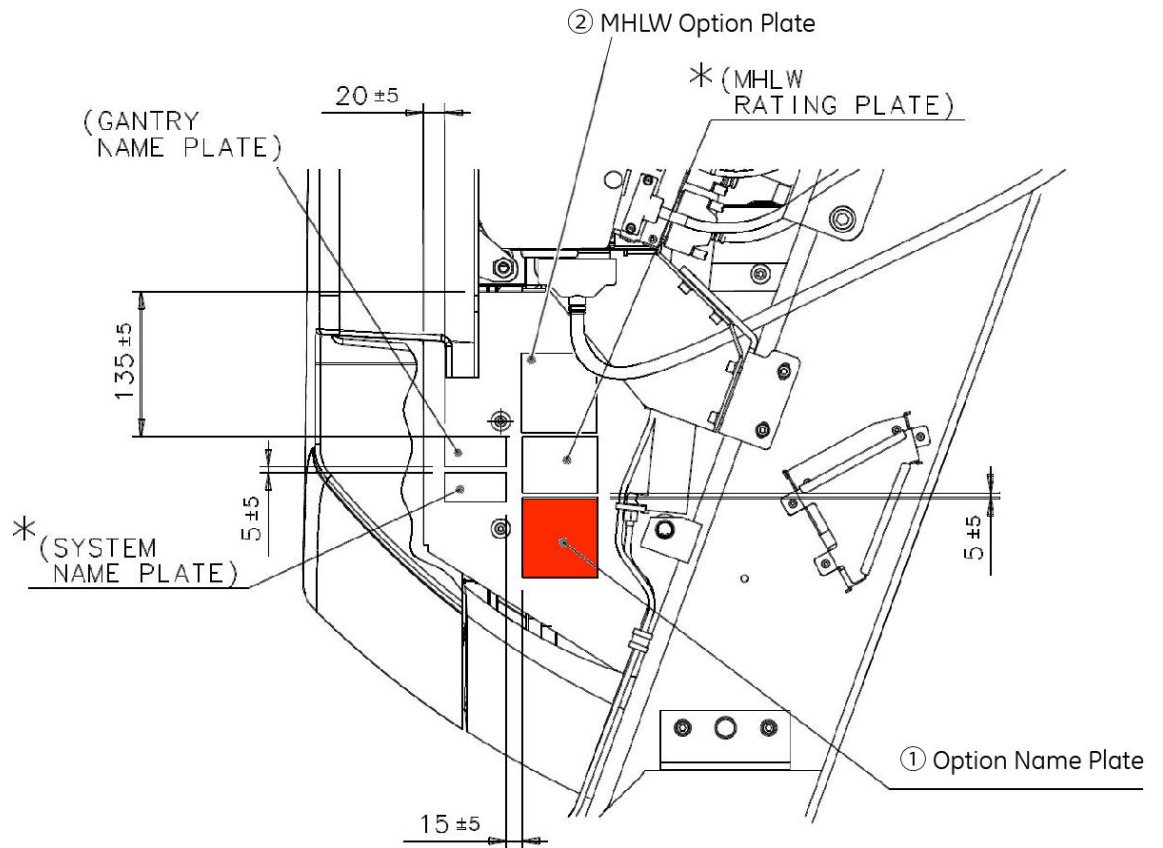
Figure 3-1 Location of Option Name Plate

- 2.) Remove the 75kVA Power Option Name Plate on the PDU rear cover. Refer to [Figure 3-2](#).

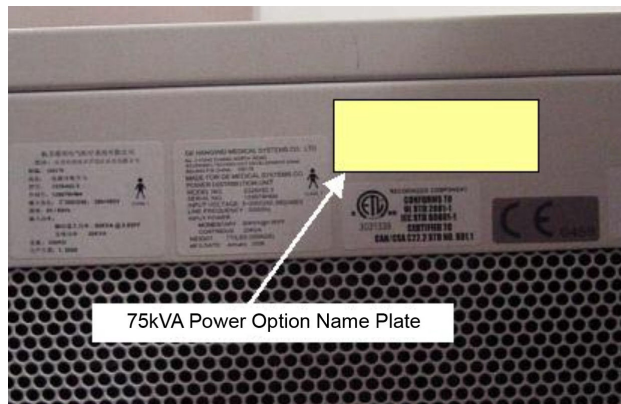
Figure 3-2 Location of Option Name Plate

## 4.2 Attach 100kVA MHLW Option Plate

- 1.) Attach the MHLW option plate near the rating plate on the right side of Gantry. Refer to [Figure 3-1](#) ②.

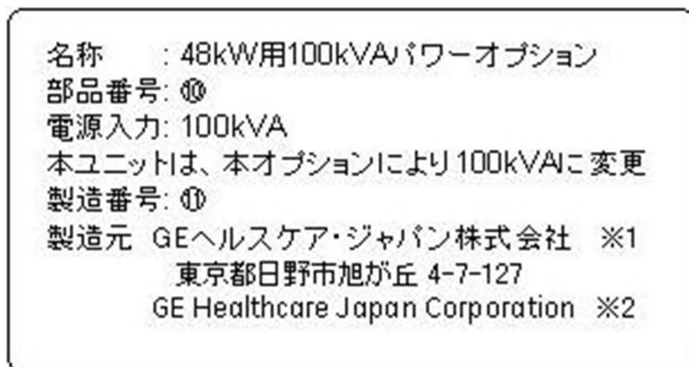
Figure 3-3 MHLW Option Plate

## 4.3 Install 100kVA Option SW Key

- 1.) Install the provided option software key.  
Go to CSD desktop -> select Configuration -> Install Options
- 2.) Reboot the system.

© 2014, General Electric Company.  
GE Medical Systems, a General Electric Company, going to market as GE Healthcare.  
3000 N. Grandview Boulevard  
Waukesha, Wisconsin 53188  
USA

[www.gehealthcare.com](http://www.gehealthcare.com)

**Imagination at work**